

# TRIGLYCERIDES

(Enzymax)

Single Reagent

## INTENDED USE

Vitro triglycerides reagent is intended for the in vitro quantitative determination of triglycerides in serum and plasma on both automated and manual systems.

## METHOD

Enzymatic colorimetric method (GPO/PAP) with glycerol phosphate oxidase and 4-aminoantipyrene.

**VITRO**  
SCIENT.

Liquid stable single reagent.


## BACKGROUND

Triglycerides are esters of the trihydric alcohol glycerol with 3 long chain fatty acids. They are the main lipids present in human plasma; the others are cholesterol, phospholipids, and non-esterified fatty acids. Triglycerides are synthesized in the intestinal mucosa by the esterification of glycerol and free fatty acids. They are then released into the mesenteric lymphatics and distributed to most tissues for storage. Triglycerides are the main storage lipids in humans, where they constitute about 95% of adipose tissue lipids. Elevated levels of triglycerides have been associated with high risk in severe atherosclerosis. High triglycerides levels and hyperlipidemia in general can be an inherited trait or can be secondary to disorders including diabetes mellitus, nephrosis, biliary obstruction, and metabolic disorders associated with endocrine disturbances<sup>1-3</sup>.

## ASSAY PRINCIPLE

Triglycerides are generally determined by a combination of hydrolysis to glycerol and free fatty acids and measurement of the amount of glycerol released. The most commonly used methods involve alkaline hydrolysis and either chemical or enzymatic measurement of glycerol. Chemical means of analysis generally rely on measurement of the product of periodate oxidation of glycerol. Eggstein and Kreutz developed an enzymatic method for measuring glycerol released from triglycerides by alkaline hydrolysis<sup>4</sup>. This method was based on the coupled reaction sequence catalyzed by glycerol kinase, pyruvate kinase, and lactate dehydrogenase. A method for complete enzymatic hydrolysis to triglycerides avoiding the need for serum pretreatment was described by Bucolo and David, using a combination of lipase and at least one proteolytic enzyme<sup>5</sup>. Wahlefeld reported that certain esterases could be combined with a lipase to achieve complete triglycerides hydrolysis<sup>6</sup>. Both methods employed a coupled enzymatic reaction sequence<sup>7</sup> to measure glycerol. Vitro triglycerides reagent combines the use of lipoprotein lipase, glycerol kinase, and glycerol phosphate oxidase with the peroxidase/4-chlorophenol/4-aminoantipyrene system of Trinder<sup>8</sup> for the measurement of triglycerides in human serum. The series of reactions involved in the assay system are as follows:

1. Triglycerides are hydrolyzed by lipoprotein lipase (LPL) to glycerol and fatty acids.
2. Glycerol is then phosphorylated to glycerol-3-phosphate by ATP in a reaction catalyzed by glycerol kinase (GK).
3. The oxidation of glycerol-3-phosphate is catalyzed by glycerol phosphate oxidase (GPO) to form dihydroxyacetone phosphate and hydrogen peroxide (H<sub>2</sub>O<sub>2</sub>).
4. In presence of peroxidase (POD), the hydrogen peroxide (H<sub>2</sub>O<sub>2</sub>) formed effects the oxidative coupling of 4 - chlorophenol and 4-aminoantipyrene (4-AAP) to form a red-colored quinoneimine dye.

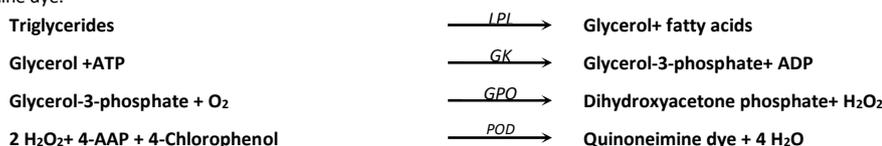

The intensity of the color produced is directly proportional to triglycerides concentration. It is determined by measuring the increase in absorbance at 500 – 550 nm.

## EXPECTED VALUES

Males:	40 – 160 mg/dl (0.45 – 1.82 mmol/l)
Females:	35 – 135 mg/dl (0.4 – 1.54 mmol/l)

For the recognition of the risk factor hyper-triglyceridemia. The following limits are recommended:

Suspicious	>150 mg/dl (1.71 mmol/l)
Elevated	>200 mg/dl (2.28 mmol/l)

Each laboratory should investigate the transferability of the expected values to its own patient population and if necessary determine its own reference range. For diagnostic purposes, the triglycerides results should always be assessed in conjunction with the patient's medical history, clinical examination, and other findings.

## REAGENTS

<b>R<sub>1</sub></b>	Triglycerides standard	200 mg/dl
	Pipes buffer, pH 7.8	50 mmol/l
<b>R<sub>2</sub></b>	p-Chlorophenol	2.0 mmol/l
	Lipoprotein lipase	1500 U/l
	Glycerolkinase	800 U/l
	Glycerol phosphate oxidase	4000 U/l
	Peroxidase	1500 U/l
	4-Aminoantipyrene	0.4 mmol/l
	ATP	0.3 mmol/l
	Mg <sup>2+</sup>	40 mmol/l
	Sodium cholate	0.2 mmol/l

### • Reagent Preparation & Stability

All reagents are ready for use and stable up to the expiry date given on label when stored at 2–8°C.

Medical Device Safety Services  
MDSS GmbH  
Burckhardtstr. 1  
30163 Hannover, Germany  
Tel.: +49 511 6262 8630  
Fax: +49 511 6262 8633

### • Serum or plasma\* from/

## SPECIMEN

- The only accepted anticoagulants are heparin and EDTA.

### Specimen Preparation & Stability For serum specimen

Patients should refrain from eating for 10 to 14 hours before blood is drawn. Samples must be drawn in a soap and glycerol free collection device. Blood should be collected by venipuncture, after the patient has been in a seated position for at least 5 minutes. Tourniquet usage should be kept to a minimum and the specimen should be allowed to clot for 30 minutes at room temperature<sup>8</sup>. The best specimen is unhemolysed serum, and should be analyzed on the day of collection. Specimens are stable for 7 days when stored at 4°C; several months at –20°C and for years at –70°C<sup>1</sup>.

## PROCEDURE

### • Manual Procedure

Wavelength	500 - 550 nm
Cuvette	1 cm light path
Temperature	20-25 or 37 °C
Zero adjustment	against reagent blank
Specimen	Serum or plasma

	Blank	Standard	Specimen
<b>R<sub>2</sub></b>	1.0 ml	1.0 ml	1.0 ml
Standard	.....	10 µl	.....
Specimen	.....	.....	10 µl

Mix, incubate for 5 minutes at 37°C or 10 minutes at 20-25°C. Measure the absorbance of specimen (A<sub>specimen</sub>) and standard (A<sub>standard</sub>) against reagent blank.

The color is stable for 60 minutes.

### • Automated Procedure

User defined parameters for different auto analyzers are available upon request

## CALCULATION

Calculate the Triglycerides concentration by using the following formulae:

Triglycerides Concentration=

Absorbance of Specimen

Absorbance of Standard  $\times$  Standard value

- **Unit conversion**  
mg/dl x 0.0114 = mmol/l

It is recommended that controls (normal and abnormal) be included in:

**QUALITY CONTROL**

- Each set of assays, or
- At least once a shift, or
- When a new bottle of reagent is used, or
- After preventive maintenance is performed or a clinical component is replaced.  
Commercially available control material with established triglycerides values may be routinely used for quality control.  
Failure to obtain the proper range of values in the assay of control material may indicate:
  - Reagent deterioration,
  - Instrument malfunction, or
  - Procedure errors.
 The following corrective actions are recommended in such situations:
  - Repeat the same controls.
  - If repeated control results are outside the limits, prepare fresh control serum and repeat the test.
  - If results on fresh control material still remain outside the limits, then repeat the test with fresh reagent.

If results are still out of control, contact Vitro Technical Services.

**INTERFERING SUBSTANCES**

- **Anticoagulants:**  
The only acceptable anticoagulants are heparin and EDTA.
- **Bilirubin:**  
No interference from free bilirubin up to level of 10 mg/dl and from conjugated bilirubin up to a level of 12 mg/dl.
- **Drugs:**  
Methyldopa and noramidopyrine causes artificially low triglycerides values at the tested drug level. For a more comprehensive review of drugs affecting triglycerides assays refer to the publication by Young<sup>9</sup>.
- **Haemoglobin:**  
No interference from haemoglobin up to a level of 600 mg/dl.
- **Lipemia:**  
Extremely lipemic samples can produce a normal triglycerides result (triglycerides greater than 3000 mg/dl).
- **Others:**  
Ascorbic acid levels higher than 2.0 mg/dl decrease the apparent triglycerides concentration significantly.

**WARNING & PRECAUTION**

- Vitro triglycerides reagent is for in vitro diagnostic use only. Normal precautions exercised in handling laboratory reagents should be followed.
- Warm up working solution to the corresponding temperature before use.
- The reagent and sample volumes may be altered proportionally to accommodate different spectrophotometer requirements.
- Valid results depend on an accurately calibrated instrument, timing, and temperature control.
- The reagent blank will not exceed an absorbance of 0.06 but don't use the reagent if it is turbid or if the absorbance is greater than 0.2 at 500 nm.
- Extremely lipemic specimens can produce a normal result. Dilute specimens 1+4 with saline and reassayed. Multiply the result by 5.

**PERFORMANCE CHARACTERISTICS**

**Imprecision**

Reproducibility was determined using in an internal protocol. The following results were obtained.

	Within Run		Between Day	
	Level I	Level II	Level I	Level II
Control				
Number of samples	40	40	40	40
Mean (mg/dl)	122	205	121	207
SD (mg/dl)	1.1	1.48	2.0	1.9
CV (%)	0.9	0.65	1.6	1.2

Manufactured in Egypt by:  
Vitro Scient

Technical Support:  
+202 26439699 [info@vitroscent.com](mailto:info@vitroscent.com)

orders:  
+202 26439698  
[order@vitroscent.com](mailto:order@vitroscent.com)

[www.vitroscent.com](http://www.vitroscent.com)

**Method Comparison**

Comparison studies were carried out using another commercially available triglyceride kit as a reference. Normal and abnormal human serum samples were assayed in

parallel and the results compared by least squares regression. The following statistics were obtained.

Number of sample pairs 45  
 Range of results 93 - 361 mg/dl  
 Slope 0.96  
 Intercept 19.5 mg/dl  
 Correlation coefficient 0.995

**Sensitivity**

The sensitivity is defined as the change of analytical response per unit change in analyte concentration at a path length of 1 cm.

When run as recommended the sensitivity of this assay is 3.0 mg/dl (0.034 mmol/l).

When run as recommended, **LINEARITY** 900 mg/dl (20.7 mmol/l).

If result exceeds 900 mg/dl (10.26 mmol/l), specimen should be diluted with 0.9% NaCl solution and reassayed. Multiply the result by the dilution factor.

**BIBLIOGRAPHY**

1. **Stein, EA (1987):** Lipids, lipoproteins, and apolipoproteins. In: Tietz NW, ed. Fundamentals of Clinical Chemistry. 3<sup>rd</sup> ed. Philadelphia: WB Saunders 448-481.
2. **Naito, HK (1984):** Disorders of lipid metabolism. In: Kaplan LA, Pesce AJ, eds. Clinical Chemistry, theory, analysis, and correlation. St. Louis: Mosby Company. 550-593.
3. **Fredrickson, DS, Levy, RI and Lees, RS (1967):** Fat transport in lipoproteins an integrated approach to mechanisms and disorders. N Engl J Med. 276: 34-43.
4. **Eggstein, M & Kreutz, F (1966):** Klin, Wschr. 44: 262-267.
5. **Bucolo, G & David, H (1973):** Clin. Chem. 19:476. **Wshahfeld, AW (1974):** Methods of Enzymatic Analysis, Bergmeyer, HU, ed pp1831, Academic press, New York.
6. **Trinder, P (1969):** Ann. Clin. Biochem. 6: 24.
7. **Henry, RJ (1974):** Clinical Chemistry: Principles and Techniques. Harper & Row, Hagerstown, MD.
8. **Young, DS (1990):** Effects of Drugs on Clinical Laboratory Tests. Third Edition. 1990: 3: 6-12.

**SYMBOL DECLARATION**

	Manufacturer
	Consult instructions for use
	Batch code (Lot #)
	Catalog number
	Temperature limitation
	In vitro diagnostic medical device
	Use by
	Caution. Consult instructions
	Keep away from light

**ORDERING INFORMATION**

REF	SIZE	REF	SIZE
13601	2 x 25 ml	13608	4 x 50 ml
13602	2 x 50 ml	13609	2 X 125 ml
13603	1 x 100 ml	13610	3 X 100 ml
13604	5 X 20 ml	13611	4 X 100 ml
13605	5 X 25 ml	13612	4 X 125 ml
13606	10 X 20 ml	13613	2 X 500 ml
13607	10 X 25 ml		

